




3nh.com

www.3nh.com


One-Button to measure



7 inch color touch screen

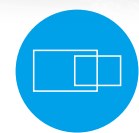

multiple measure model

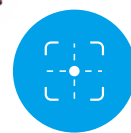

camera locating

YS6020-Benchtop Grating Spectrophotometer

The spectrophotometer YS6020 is also widely used in the control of color difference quality of various products. YS6020 is equipped with color management software, which can be connected to a computer to achieve more functional expansion. The instrument has stable performance, accurate color measurement and powerful functions. It is widely used in scientific research institutions and laboratories in industries such as plastic electronics, paints and coatings, textile printing and dyeing, printing paper products, automotive, medical, cosmetics and food.

APPLICATION INDUSTRY

YS6020 benchtop spectrophotometer is designed for accurate analysis and transmission of laboratory color and appearance: It's widely used in plastics, paints, inks, textile and garment printing, dyeing, printing paper, automotive, medical, cosmetics, food industries, scientific research institutions and laboratories.



Automobile



Leather



Plastics




Paint



Food stuff



Laboratory

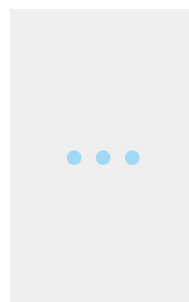

Others

PRODUCT FEATURES

- 1、 High configuration of hardware: 7 inches TFT Color Capacitive Touch-screen Display; Bluetooth; Concave Grating.
- 2、 Double Array 256 Image Element CMOS Sensor; Long life-span stable xenon lamp.
- 3、 With reflective and transmissive spectrum, accurate Lab value, good to calculate color formula and do precise color transmission.
- 4、 Auto identify measuring aperture. Freely switchable between 3 measuring apertures: Φ 25.4mm/8mm/4mm. Users also can customize apertures.
- 5、 Built-in temperature sensor to monitor and compensate the measuring temperature to ensure the measurement more precision.
- 6、 Wavelength range 360nm – 780nm. Built-in 400nm cut off, more professional in UV measurement.
- 7、 Independent light source detector, continuously monitor the condition of light sources to ensure the light source reliable.
- 8、 Multiple measurement modes: Quality Management Mode, Sample Mode; Meet more users' requirement.
- 9、 A variety of optional accessories: Reflection sample holding tool, transmission fixture, micro 4mm aperture transmission test components, instrument inversion test components, applicable to more working conditions.
- 10、 Big capacity data storage, for 20000 pieces test result.
- 11、 Built-in camera locating.
- 12、 More powerful extended functions at the PC software.




PRODUCT FEATURES



plastic sample

Reflect: di:8°, de:8°(diffused illumination, 8-degree viewing angle) Conforms to ISO7724/ 1、 CIE、 ASTM、 DIN and JIS

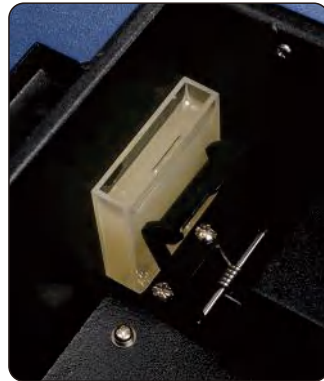

Transmissive samples

Transmittance: di:0°, de:0° (diffuse illumination: 0° viewing); Conforms to ISO、 CIE、 ASTM 、 DIN.



Liquid, Powder, Solid samples



Measuring Aperture

Auto identify measuring aperture. Freely switchable between 3 measuring apertures: Φ 25.4mm/8mm/4mm.

OPTIONAL ACCESSORIES

The instrument is equipped with a variety of optional accessories: micro printer, transmission test fixture assembly, instrument inversion test fixture, petri dish; Meet customer multi-sample measurement.



Microprinter



Instrument upside down test fixture




Transmission test fixture assembly







Petri dish

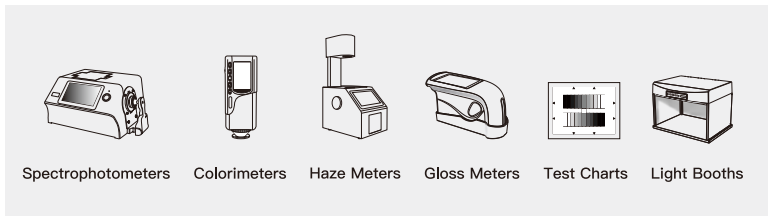Model: YS6020(Pulsed xenon lamp)
Optical Geometry: Reflect: di:8°, de:8°(diffused illumination, 8-degree viewing angle);SCI (specular component included)/SCE (specular component excluded);Include UV / excluded UV light source; Transmittance: di:0°, de:0° (diffuse illumination: 0° viewing);SCI (specular component included)/SCE (specular component excluded);Include UV / excluded UV light source; Haze(ASTM D1003); Conforms to CIE No.15, GB/T 3978,GB 2893,GB/T 18833, ISO7724/1, ASTM E1164, DIN5033 Teil7
Integrating Sphere Size: Φ 154mm
Light Source Device: 360nm-780nm Xenon Lamp, 400nm cut-off Xenon Lamp
Spectrophotometric Mode: Concave Grating
Sensor: 256 Image Element Double Array CMOS Image Sensor
Wavelength Range: 360-780nm
Wavelength Pitch: 10nm
Semiband Width: 5nm
Reflectance Range: 0~200%
Measuring Aperture: Reflective: Φ30mm/Φ25.4mm, Φ10mm/Φ8mm, Φ6mm/Φ4mm;Transmissive: Φ30mm/Φ25.4mm; Remark: 1. Automatic identification of switch caliber 2. Customized Configuration caliber and lens position
Specular Component: Reflectance: SCI&SCE Transmittance: SCI&SCE
Color Space: CIE Lab, XYZ, Yxy, LCh, CIE LUV, Hunter LAB, Munsell, s-RGB, HunterLab, DIN, βxy
Color Difference Formula: ΔE*ab, ΔE*uv, ΔE*94, ΔE*cmc(2:1), ΔE*cmc(1:1), ΔE*00v, ΔE(Hunter), DIN ΔE99
Colorimetric Index: WI (ASTM E313, CIE/ISO, AATCC, Hunter), YI (ASTM D1925, ASTM 313), MI (Metamerism Index), 8° glossiness Staining Fastness, Color Fastness, Color Strength, Opacity,Gardner Index, APHA/Hazen/Pt-Co Index,555 Index, Haze(ASTM D1003) Pharmacopoeia (China, Europe, the United States), EBC (Beer Color), ICUMSA Color (Sugar Color)
Observer Angle: 2° & 10°
Illuminants: D65, A, C, UV, D50 , D55, D75, F1, F2, F3, F4, F5, F6, F7, F8, F9, F10, F11, F12 , CWF,DLF,TL83,TL84,TPL5,U30
Displayed Data: Spectrogram/Values, Chromaticity Values, Color Difference Values/Graph, Pass/Fail Result, Color Offset
Measuring Time: About 2.4s (Measure SCI & SCE about 5s)
Repeatability: Spectral reflectance: Φ25.4mm/SCI, Standard deviation within 0.07% (400 nm to 700 nm: within 0.06%)Chromaticity value: Φ25.4mm/SCI, Standard deviation within ΔE*ab 0.015 (When a white calibration plate is measured 30 times at 5 second intervals after white calibration)Spectral Transmittance: Φ25.4mm/SCI, Standard deviation within 0.07% (400 nm to 700 nm: within 0.07%)Chromaticity value: Φ25.4mm/SCI, Standard deviation within ΔE*ab 0.018 (When a white calibration plate is measured 30 times at 5 second intervals after white calibration)
Inter-instrument Error: Φ25.4mm/SCI, Within ΔE*ab 0.15 (Average for 12 BCRA Series II color tiles)
Size: 370×300×200mm
Weight: About 9.6kg
Power Supply: DC 24V, 3A Power Adapter
Light Source Device Life: 5 years, more than 3 million times measurements.
Display: 7" TFT Capacitive Screen-touch Display
Data Port: USB & Print serial port
Data Storage: Standard 2000 Pcs, Sample 20000 Pcs
Network calibration: Supports Netmetric
Language: Simplified Chinese, Traditional Chinese, English , Russian, Italian
Standard Accessory: Power Adapter, User Guide, CD Disk(PC Software), USB cable, Standard Calibration Board, Black Calibration Cavity, Transmission black baffle, Sample holder, 25.4 caliber, 8 caliber, 4 caliber
Optional Accessory: Micro-printer, Transmissive Test Clamp Component, Micro Aperture(4mm) transmission test clamp component, Instrument inversion components ,Culture dish
Notes: The specifications are subject to change without notice.

GUANGDONG THREENH TECHNOLOGY CO., LTD.



★ CONTACT US

- 
[web:www.3nh.com](http://www.3nh.com)

[Email:3nh@3nh.com](mailto:3nh@3nh.com)
- 
[Tel:0086-020-82880288](tel:0086-020-82880288)
- 
 Add: 6-8th floors, Building B33, Low Carbon Headquarters Park, Xincheng Road No.400, Zengcheng District, Guangzhou, Guangdong Province, China